

24-Hour Ambulatory Blood Pressure Monitor

SunTech Medical
Brilliant Blood Pressure Solutions™

BravoTM **24-HR ABP**



The best value in ambulatory technology is also the easiest to use.

- Fast setup
- Low power consumption
- Compact and lightweight



Bravo™ - Ambulatory Blood Pressure Monitoring System



The New Standard in ABP

www.SunTechMed.com
sales@SunTechMed.com


The *Bravo* ambulatory blood pressure monitor has a compact and efficient design providing a cost-effective solution for monitoring your patients' blood pressures outside the clinical environment. The *Bravo* monitor interfaces with AccuWin Pro™ PC software for quick programming, data retrieval and report generation.

FEATURES AND BENEFITS

Affordable:	Exceptional value conserves limited budget resources.	Data Security:	BP data saved on flash memory means no lost data and no repeated studies.
Quick Setup:	Oscillometric technology makes patient set-up easy, saving your staff time.	Low Power Consumption:	Requires only two AA batteries, making the <i>Bravo</i> much less expensive to operate than most alternatives.
Ready to Use:	Use the system right away with included ABP cuff, accessories and batteries.	Patient Comfort:	Adjustable maximum and dynamic inflate pressures enhance patient comfort and compliance.

ACCUWIN PRO

- Create a monitoring schedule tailored to an individual patient with up to 3 specific time periods.
- Generate customized clinical reports with clearly displayed ABP data in numerous graphical and/or tabular formats.
- Reported diagnostic information includes all collected BP readings, heart rate, mean arterial pressure, pulse pressure, BP load and averages calculated for day, night and overall study.

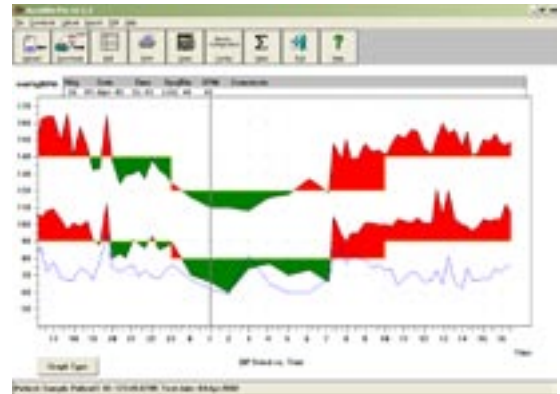

WHAT'S INCLUDED

- *Bravo* ABP Monitor
- *AccuWin Pro* Software CD
- Monitor pouch, waist belt and shoulder strap
- 4 AA batteries
- SunTech Adult ABP cuff (27-42 cm range)
- Patient Diaries
- User Manuals & Quick Guide

SPECIFICATIONS

Technique:	Oscillometry with step deflation	Event Key:	Programmable setting allows patient to initiate additional readings.
Validations:	ANSI/AAMI SP10:2002	Power Source:	3 VDC - Two 1.5V "AA" Alkaline or NiMH rechargeable batteries
BP Range:	25-260 mmHg	Dimensions:	Approx: 12 X 7 X 3 cm
HR Range:	40-200 bpm	Weight:	265 grams, including batteries
Data Memory:	Flash memory stores up to 250 readings	PC Requirements:	9 pin serial port (USB adapter optional) Windows 98 and later.
Sampling Periods:	3 independently programmable time periods		
Options:	Language versions available		
Warranty:	1 year standard warranty		



507 Airport Blvd., # 117
Morrisville, North Carolina 27560-8200 USA
Phone: 919.654.2300
Fax: 919.654.2301
E-mail: sales@SunTechMed.com

Oakfield Industrial Estate
Stanton Harcourt Road
Eynsham, Oxfordshire OX29 4TS England
Phone: +44 (0) 1865 884 234
Fax: +44 (0) 1865 884 235

82-0024-00 Rev. D